

# AMRAD™

Overhead Tube Support

EXAMINE  
WHAT AMRAD™  
CAN DO FOR YOU.

Take a closer look

# AMRAD ADVANTAGE OVERHEAD TUBE SUPPORT.

## Reliable

AMRAD™ is known for innovative development and quality workmanship in radiographic systems and components. AMRAD™ products help physicians provide top-quality health care for their patients every day. Each component has been designed for long-term durability under heavy patient volume conditions. Backed by a five-year parts warranty, AmradAdvantage radiographic equipment offers exceptional performance, reliability and long-term investment value.

## Versatile

The Overhead Tube Support (OTS) is ideal for sophisticated imaging environments. Gliding smoothly through a complete range of motion, the OTS offers flexibility and efficiency for performing a full array of anatomical studies with precision and ease.

Longitudinal rails afford virtually unlimited travel for complete tabletop coverage. The extensive vertical travel permits exams ranging from skull through weight bearing. The tube head rotates 300 degrees around the telescopic vertical column and 240 degrees around the horizontal axis with positive detents every 90 degrees, simplifying table, off-table and upright procedures. Tube trunnion rings permit unrestricted cross-table imaging.

## Accurate

The fully counterbalanced, telescopic column with double cable safety system allows zero drift, ensuring safe and accurate operation.


The overhead structure rides on stainless steel precision bearings and hardened steel tracks to ensure sturdy, vibration-free positioning. Electromagnetic locks control all movements. Individual lock releases are conveniently located on the user-friendly operator handle including a center mounted, single-button all-lock release for quick positioning of the tube in emergency situations. Solid-state liquid crystal displays on the operator handle provide clear visual indications of the Source Image Distance to the tabletop and bucky, as well as the tube angle. Designed to accommodate high patient volumes, the OTS maximizes efficiency and simplifies precise tube positioning.



Upright procedure shown with OTS and self-centering wall cassette stand.



Weight bearing procedures shown with OTS and self-centering wall cassette stand.

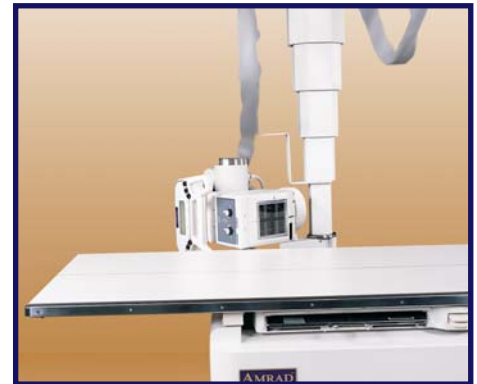

Cross-table procedures shown with OTS and elevator table.

## Worldwide Support

AMRAD™ radiographic products are backed by an exclusive, worldwide dealer network for support and service. This exceptional group of professionals is dedicated to providing you with unequalled products and services.